



Introducing

G3 X

The First Standard Size
Foam-Free PAP
Available in the U.S.



Five Year Warranty

Long-term support and reliability for your patient's entire therapy journey.

Foam Free

The first standard-size PAP available in the U.S. engineered without polyurethane sound-abatement foam while maintaining quiet operation.

Proven Algorithm

Demonstrated accurate identification of obstructive events through independent bench study.



Smart Features

Custom pressure response and Smart A continuously responds to each patient's unique needs for unparalleled comfort, delivering precision care.

Seamless Care Management

The integrated cloud platform enables direct provider-patient communication and actionable insights that transform care delivery.



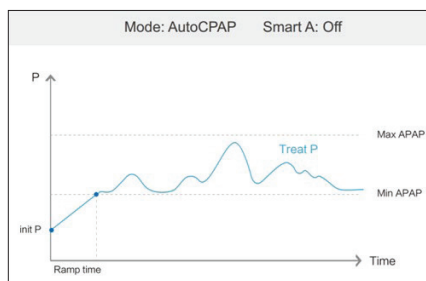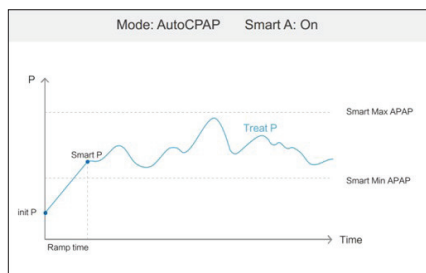

Smart A Custom Pressure Response for Precision Care

A patient's condition is always changing. Smart A fine-tunes the Auto Algorithm to respond to the patient's pressure needs to support comfort and provide success with PAP therapy.

The Smart A algorithm uses data from a rolling 5-effective treatment sessions. An effective treatment session is defined as at least 4 hours of use and a controlled leak (<90 lpm). The algorithm looks at the 50th and 95th percentile to determine where Auto min and Auto max should be adjusted.

How G3 X Stacks Up

	G3 X	G3	AirSense™* S11	AirSense™* S10
Brand	React Health	React Health	ResMed	ResMed
Warranty	5 years	2 years	2 years	2 years
Foam Free	Yes	No	No	No
Precision care with custom pressure response	Yes	No	No	No
Connectivity	4G, QR Code, SD Card	4G, QR Code, SD Card	4G, SD Card	4G, SD Card
Connected Platform	Yes	Yes	Yes	Yes
Patient App	Yes	Yes	Yes	Yes
Directly communicate to patients through the patient app	Yes	Yes	No	No

REPLACEMENT PARTS

G4600	G3 X APAP Device with Cellular Modem-HT
G46002	G3 X APAP Device with Cellular Module without HT
LG34510	G3 Replacement Water Chamber
LG3HT	G3 Integrated Heated Tubing
UN100	Ultra Noir CPAP Standard Tubing (Non-Heated) CS of 50

LG3F002	G3 Disposable Ultra Fine Filter
LG38008	G3 Reusable Filter - Box of 12 Filters (Black)
LG3B	G3 Repl Carry Bag
LPG3	G3 Power Supply Adapter
LPAC	G2/G3 Power Cord
PRTG301	G3 Repl Filter Cover- Black

ReactHealth.com

© 2025, 3B Medical, Inc. dba React Health, all rights reserved.

Caution: U.S. Federal law restricts this device to sale by or on the order of a physician. Technical specifications may change.

*AirSense™ is a registered trademark of ResMed family of companies.



MI.MK.0413 Rev A