

# Ventec Clinical Support

Every VOCSN comes with personalized support from the Ventec Team

## Setup

We partner with your team to set up VOCSN and help you knowledgeably transition patients.

## Train

We personally train every member of your team to use each VOCSN therapy.

## Support

Ventec customer support with experienced doctors, respiratory therapists, and engineers is just a phone call away.



# VOCSN

Five therapies. One device.



**Customizable:** Get all therapies or just the mix of therapies needed

Configuration	Ventilation	O2 Concentrator	Cough	Suction	Nebulizer	High Flow Nasal Cannula Therapy	High Pressure External O2 and FiO2 Monitor	Low Pressure External Oxygen
<b>V+O+C+S+N+Pro</b> (PRT-00853-000)	✓	✓	✓	✓	✓	✓	✓	✓
<b>V+O+C+S+N</b> (PRT-01126-000)	✓	✓	✓	✓	✓	✓	—	✓
<b>V+C+S+N+Pro</b> (PRT-01172-000)	✓	—	✓	✓	✓	✓	✓	✓
<b>V+C+Pro</b> (PRT-00973-000)	✓	—	✓	—	—	✓	✓	✓
<b>V+C</b> (PRT-01125-000)	✓	—	✓	—	—	✓	—	✓
<b>V+Pro</b> (PRT-01198-000)	✓	—	—	—	—	✓	✓	✓



VENTILATION



OXYGEN



COUGH



SUCTION



NEBULIZER

Integration makes care quick and easy

## HOSPITAL TO HOME

VOCSN



VENTILATION

### Critical Care Ventilator

VOCSN is one of the few portable ventilators to achieve the rigorous critical care standard for safety and accuracy.\* The ventilator provides invasive, noninvasive, and mouthpiece ventilation and delivers a comprehensive set of modes and settings to meet patient needs.



OXYGEN

### 6 L/min Equivalent Internal Oxygen Concentrator

The VOCSN internal oxygen concentrator and Oxygen Direct® system deliver the equivalent of 6 L/min of oxygen or up to 40% FiO2. External high pressure and low pressure oxygen sources can be connected when needed.



COUGH

### Touch Button Cough®

Touch Button Cough therapy is activated in seconds rather than minutes, without changing the circuit. Breath Sync monitors patient breathing and triggers a cough at a natural point in the breathing cycle.



SUCTION

### Hospital Grade Suction

Hospital grade suction system provides consistent high flows throughout the entire suction experience. VOCSN suction therapy is up to three times quieter than traditional portable suction machines.



NEBULIZER

### High Performance Nebulizer

VOCSN automatically compensates for the airflow from the nebulizer drive to ensure accurate ventilation and comfortable breathing. VOCSN records each medication use and turns off the nebulizer when the therapy is complete.

**VENTILATOR**  
12.4 lb

**COUGH ASSIST**  
9.4 lb

**SUCTION**  
11.4 lb

**OXYGEN CONCENTRATOR**  
18.4 lb

**NEBULIZER**  
3.3 lb

VS



**55lb** 5 MACHINES  
MULTIPLE CIRCUITS

**12.5 - 18lb** 1 MACHINE  
1 CIRCUIT

Depending on configuration

\*ISO 80601-2-12



## Benefits of Integration



### SIMPLE

- Ventec One-Circuit®
- Intuitive touchscreen operation
- Quick and easy setup

**1** device to learn  
device to clean  
alarm to silence



### MOBILE

- Five therapies in one device
- 70% lighter and smaller than five devices
- Up to 9 hours of on-board battery


**1** device to transport  
device to charge  
power supply to carry



### CARE CHANGING

- Touch Button Cough
- Uninterrupted ventilation between therapies
- Designed to reduce risk of circuit misconnects

**1** circuit for all therapies  
touch to operate  
device to monitor

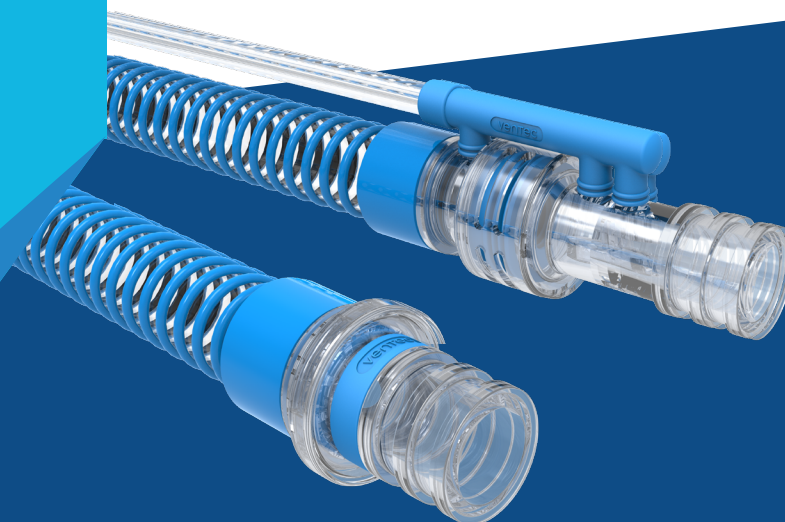

## Ventec One-Circuit

Ventec One-Circuit allows switching between therapies with the touch of a button. The redesigned valves make Touch Button Cough therapy possible and the integrated oxygen tube enables Oxygen Direct. Rather than several cumbersome tubes connecting to separate devices, the Ventec One-Circuit eliminates clutter to create one easy-to-manage system for peace of mind.