
REACTHEALTH

## LUNA G3 BILEVEL ST (LG3700 / LG3800)



The Luna G3 Bilevel ST 30Vt model is a Bilevel PAP device intended to provide non-invasive ventilation in the home or hospital/institutional environment

With its Target Tidal Volume function, excellent comfort features and effective performance, the Luna G3 Bilevel ST 30Vt offers each patient personalized ventilation support

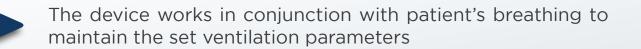

The device detects leaks in real time and adjusts the baseline to ensure correct triggering and cycling





Broad Range of Alarm Settings



Leak Compensation



Target Tidal Volume Function



Efficiencies in Therapy with S/T Mode



Inspiratory / Expiratory Sensitivity

## Luna G3 Specifications Overview

Specification	Luna G3 Bilevel ST 30 Vt (LG3800)	Luna G3 Bilevel 25 A (LG3700)
Dimensions (L x W x H)	10.4" x 5.7" x 4.5"	10.4" x 5.7" x 4.5"
Weight	3.75 lbs	3.75 lbs
Humidification	Standard, Integrated	Standard, Integrated
Heated Tubing	Standard	Standard
Connectivity	iCode, Integrated Cellular Modem*, QR Code, SD Card, Wi-Fi*	iCode, Integrated Cellular Modem*, QR Code, SD Card, Wi-Fi*
SD Data Capacity	Minimum 12 months	Minimum 12 months
Filters	Reusable or Disposable	Reusable or Disposable
Mode of Operation	CPAP, S, T, S/T	CPAP, AutoCPAP, S, AutoS
Ramp Time	0 to 60 min (5 min increments)	0 to 60 min (5 min increments)
Inspiration Time	0.1 to 3 seconds (0.1 sec increments)	
Ventilation Pressure	4 to 30 cmH2O (0.5 cmH2O increments) EPAP: 4 to 25 cmH2O (0.5 cmH2O increments) CPAP max pressure 20 cmH2O	4 to 25 cmH2O
Weight Restrictions	Over 66 lbs	Over 66 lbs
Sound Pressure Level dBA	< 26 dB, when device is working at the pressure	< 26 dB, when device is working at the pressure
Sound Power Level dBA	< 34 dB, when the device is working at the pressure	< 34 dB, when the device is working at the pressure

\*Additional fees are required to support cellular and/or Wi-Fi functionality

MI.MK.0337.REV.A.



change.